

Made for
iPhone | iPad | iPod

Works with
android



User guide

Beltone Custom hearing aids



Beltone **Rely™**

Hearing aid information

Left hearing aid		Right hearing aid	
Serial number		Serial number	
Model number		Model number	
Battery type	☐ CIC size 10A ☐ ITC size 312, ITE size 312 ☐ ITE size 13		

Program	Beep	Description
1	One beep	
2	Two beeps	
3	Three beeps	
4	Four beeps	

Table of Contents

Introduction	5
Your hearing aid	7
How to get your hearing aid ready for use	9
How to place the hearing aids in your ears	14
How to remove the hearing aids from your ears	15
How to use your hearing aids	16
Using a telephone	20
Advanced options	24
Wireless accessories	28
How to clean and maintain your hearing aids	30
General warnings and precautions	33
Hearing aid expectations	35
Troubleshooting	36
Tinnitus Management	37
Regulatory information	45
Additional information	50

Introduction

Thank you for choosing Beltone hearing aids. We recommend that you use your hearing aids every day. This way you will fully benefit from them.

NOTE: Read this booklet carefully before you start using your hearing aids.

Intended use

Generic air-conduction hearing aids are wearable sound-amplifying devices intended to compensate for impaired hearing. The fundamental operating principle of hearing aids is to receive, amplify, and transfer sound to the eardrum of a hearing impaired person.

Contra-indications

A hearing care professional should advise a prospective hearing aid user to consult promptly with a licensed physician (preferably an ear specialist) before dispensing a hearing aid, if the hearing aid dispenser determines through inquiry, actual observation, or review of any other available information concerning the prospective user, that the prospective user has any of the following conditions:

- Visible congenital or traumatic deformity of the ear.
- History of active drainage from the ear within the previous 90 days.
- History of sudden or rapidly progressive hearing loss within the previous 90 days.
- Acute or chronic dizziness.
- Unilateral hearing loss of sudden or recent onset within the previous 90 days.
- Audiometric air-bone gap equal to or greater than 15 dB at 500 Hertz (Hz), 1000 Hz, and 2000 Hz.

- Visible evidence of significant cerumen accumulation or a foreign body in the ear canal.
- Pain or discomfort in the ear.

Side effects

If you experience side effects, contact your hearing care professional. Possible side effects from wearing a hearing aid may be:

- Dizziness
- Tinnitus
- Perceived worsening of hearing loss
- Nausea
- Skin reaction
- Ear wax accumulation



For devices including a Tinnitus Breaker Pro module

The Tinnitus Breaker Pro module is a tool to generate sounds to be used in a Tinnitus Management Program to temporarily relieve patients suffering from Tinnitus.

Your hearing aid


Your hearing aid – Completely-In-the-Canal (CIC)

1. Push button (optional)
2. Battery door and On/Off switch
3. Removal cord
4. Sound outlet
5. Wax filter
6. Vent
7. Microphone inlets
8. Manufacturer
9. Model designation and serial number
10. Antenna

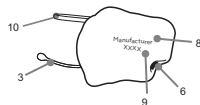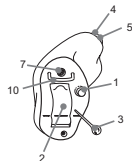

Your hearing aid – In-The-Canal (ITC)/In-The-Ear (ITE)


1. Microphone sound inlets
2. Push button (optional)
3. Battery door
4. Volume control (optional)
5. Vent

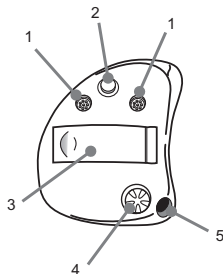


How to get your hearing aid ready for use

Replacing the battery

1. Open the battery door completely by using your fingernail. Remove the used battery if present.

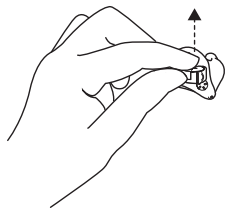

2. Prepare the new battery. Remove the protective foil to activate the battery. Wait for 2 minutes before inserting the battery into the hearing aid.

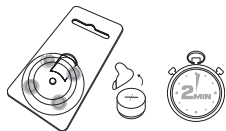

3. Insert the new battery with the positive side in the correct position. Always insert the battery in the door: never directly into the hearing aid. Gently close the battery door.



When the hearing aid is not in use, open the battery door to turn it off. Open the battery door completely to allow moisture to evaporate and prolong the hearing aids' life span.

If the hearing aids are experiencing frequent loss of connection to wireless accessories, contact your hearing care professional for a list of low impedance batteries.

Keep spare batteries on hand.



Battery warnings

- Never put a hearing aid battery in your mouth.
- Keep hearing aid batteries away from pets, children, and people with cognitive, intellectual or mental health challenges.
- Never swallow batteries nor place them inside any part of the body, as the battery can cause serious injuries. If a battery has been swallowed or placed inside any part of the body, seek immediate medical attention.
- If the batteries are not inserted correctly, the device will not work and the batteries may build up heat. If this happens, please remove the batteries.
- Batteries are harmful for the environment. Therefore, never try to burn them. Dispose of your used rechargeable hearing aids according to your country's regulations or return them to your hearing care professional.
- Batteries may leak. If you are not going to use your hearing aids for a few days, you **MUST** remove the batteries.
- Battery leakage can cause chemical burns. If you are exposed to a battery leakage, rinse immediately with lukewarm water. If you get chemical burns, redness or skin irritation, seek immediate medical attention.
- Never try to charge zinc-air batteries as this can cause leakage or a small explosion.

NOTE: Always use new zinc-air batteries that have a minimum remaining shelf life of one year.

Low battery warning

When the batteries are low on power, your hearing aids reduce the volume, and play a melody every 15 minutes, until they are empty and turn off.

NOTE: Keep spare batteries on hand.

Low battery alert when paired with wireless accessories (optional)

NOTE:

- Your batteries will drain faster when you use wireless functions like streaming from your smart device or from your TV with our TV streamer. As the battery power declines, the wireless functions stop working. A short melody will play every five minutes to let you know that the battery power is low. The table below shows the functionality with different battery charge levels.
- If the hearing aids are experiencing frequent loss of connection to wireless accessories, contact your hearing care professional for a list of low impedance batteries.

Battery level	Signal	Hearing aid	Remote control	Streaming
Fully charged		✓	✓	✓
Low	 4 even tones	✓	✓	x
Depleted	 3 even tones and 1 longer tone	✓	x	x

These will work again when you insert a new battery.

How to place the hearing aids in your ears


How to tell left from right

Your hearing aids are colour-coded. Left = blue. Right = red. If your hearing aids are not colour coded, ask your hearing care professional to add colour coding.

Inserting the hearing aids

Insert Completely-In-The-Canal (CIC), In-The-Canal (ITC), and In-The-Ear (ITE)

1. Hold the hearing aid between your thumb and index finger, either above and below or on the sides.
2. Place the sound outlet portion into your ear canal. Turn the top part of the hearing aid gently backwards and forwards so that it tucks behind the fold of skin above your ear canal.
3. Insert the hearing aid into your ear canal. Opening and closing your mouth may ease insertion.

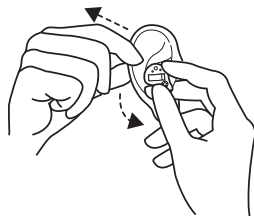


By experimenting, you may discover an easier method. With proper insertion, hearing aids should fit snugly but comfortably.

NOTE: It may be helpful to pull your ear up and outward with your opposite hand during insertion.

How to remove the hearing aids from your ears

Remove Completely-In-Canal hearing aids

- Hold the edges of the hearing aid with your thumb and index finger. Pull out and slightly upward, while slightly rotating your hand forward.

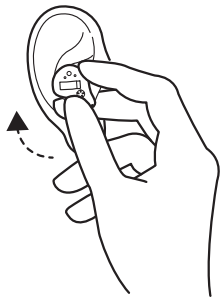

Remove In-The-Canal (ITC) and In-The-Ear (ITE)

Hold the edges of the hearing aid with your thumb and index finger. Pull out and slightly upward, while slightly rotating your hand forward.




How to use your hearing aids

Operation of the hearing aid

The volume control (optional) allows you to adjust the volume of your hearing aids to your liking.

1. To increase the volume, turn the volume control up (+).
2. To decrease the volume, turn the volume control down (-).

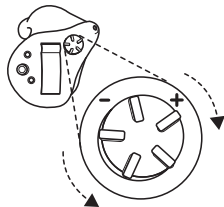

When you change the volume, the hearing aid responds with a beep. When you reach the upper or lower limits, the hearing aid responds with a low-pitch beep.

NOTE: If you have two hearing aids with the Synchronised Volume Control function enabled, volume control adjustments to one hearing aid automatically repeats in the second hearing aid. When you change volume on one of the hearing aids, it responds with one or more beeps. A beep in the second hearing aid follows.

If your hearing aids support wireless remote control, you may use the Beltone Direct Remote Control 2 or the Beltone HearMax™ app to adjust e.g. the volume.

Push button (optional)

Your hearing aid has a push button allowing you to use up to four different listening programs. The list on page 2 tells which programs have been enabled.

1. Push the program button to switch between programs.
2. You will then hear one or more beeps. The number of beeps indicates which program you have selected (one beep = program one, two beeps = program two, etc.).
3. When you turn the hearing aids off and then back on, they always return to the default setting (program 1 and preset volume).

NOTE: If you have two hearing aids with the Synchronised Push Button enabled, program changes to one hearing aid automatically repeats in the second hearing aid. The same number of confirmation beeps will follow in the second hearing aid. This Synchronised Push Button can also be configured to allow one side to control volume increase and the other to control volume decrease. The volume changes to one hearing aid are repeated on the other side to keep levels the same.

Telecoil

(Optional for ITC, ITE models only)

Your hearing aid may have a telecoil. The Telecoil function may help to improve understanding of speech with Hearing Aid Compatible (HAC) telephones and in theatres, cinemas, houses of worship, etc. that have a teleloop installed.

When you select the Telecoil program, your hearing aid picks up signals from the teleloop or HAC telephone. Your hearing care professional can activate the Telecoil program.

NOTE:

- The telecoil does not work without a teleloop (that is, an induction loop) or a HAC telephone.
- If you are having trouble hearing with the telecoil, ask your hearing care professional to adjust the program.
- If there is no sound from your hearing aids in a teleloop system with an active Telecoil function, the teleloop system may not be turned on or may not be operating correctly.
- The sound from the teleloop and the hearing aids' microphones can be mixed according to your preference. Ask your hearing care professional for more details.

How to use the telecoil program

To use teleloop systems, follow these steps:

1. Switch your hearing aid to the Telecoil program.
2. Find a good spot. Reception is not clear in all locations, it depends on the teleloop. Look for signs or find another spot.
3. If needed, adjust the volume.
4. When you leave, switch to program 1.

HAC telephone

Some smartphones are hearing aid compatible (HAC). The HAC phone establishes a small hearing loop that your hearing aids can connect to. The telecoil picks up the HAC phone's signal and converts it to sound.

To use a HAC phone, follow these steps:

1. Switch your hearing aid to the Telecoil program.
2. Pick up the phone and make a call or answer a call.
3. Hold the phone close to the hearing aid and tilt it slightly outwards.
4. Listen to the dial tone and move the telephone to get the best reception.
5. If needed, adjust the volume.
6. When you hang up, switch back to your preferred program.

NOTE:

- If the phone has a poor telecoil signal, use the microphone program. To avoid whistling, do not hold the handset too tightly against your ear.
- Ask your hearing care professional to enable the Telecoil programme in your hearing aids.
- If you see a "M3", "M4", "T3", or "T4" on the box, then the smartphone is HAC compliant. If you find it difficult to obtain a good result while using your smartphone, your hearing care professional will be able to give you advice on available wireless accessories to enhance listening capabilities. Ask your hearing care professional for advice regarding HAC smartphones.

Using a telephone

Your hearing aid allows you to use your telephone as you normally do. Finding the optimal position for holding the phone may require practice.

The following suggestions may be helpful:

1. Hold the telephone up to your ear canal as usual.
2. If you hear whistling, try holding the telephone in the same position for a few seconds. The hearing aid may be able to cancel the whistling.
3. You can also try holding the telephone slightly away from the ear.

NOTE: Depending on your needs, your hearing care professional may activate a function specifically for telephone use.

Mobile phones

Your hearing aids comply with the most stringent Standards of International Electromagnetic Compatibility. Any degree of disturbance can be due to the nature of your particular mobile phone or of your wireless telephone service provider.

NOTE:

- If you find it difficult to get a good result while using your mobile phone, your hearing care professional can give you advice on available wireless accessories to enhance listening capabilities. See also "Wireless accessories", page 28.

Auto-Phone (optional)

If you place a magnet on the telephone receiver, your hearing aids will automatically switch the telephone program on when the receiver is close to your ear. When you remove the receiver from your ear, the hearing aids automatically return to the previous listening program.

NOTE: Ask your hearing care professional to enable Auto-Phone as one of your programs.



Auto-Phone warnings

- If a magnet is swallowed, seek immediate advice from a medical practitioner.
- Keep magnets out of reach of pets, children and mentally disabled persons.
- The Auto-Phone magnet may affect sensitive medical devices/electronic systems. Seek advice from the manufacturers regarding appropriate safety measures when using the Auto-Phone solution near the sensitive

device/equipment (pacemakers and defibrillators) in question. If the manufacturer cannot issue a statement, we recommend keeping the magnet or a telephone equipped with the magnet 30 cm (12") away from magnetically sensitive devices (e.g. pacemakers).


Auto-Phone precautions

- If you experience frequent signal loss or noise during calls, move the magnet to another place on the telephone receiver.
- Only use magnets supplied by Beltone.


Place the Auto-Phone magnet

Place the magnet on your telephone receiver to allow operation of the Auto-Phone function. In order to place the magnet properly:


1. Clean the telephone thoroughly

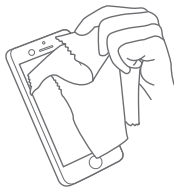

2. Remove foil from magnet

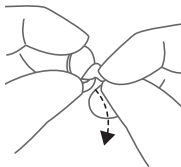

3. Place the magnet

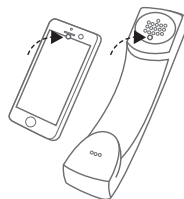

NOTE: If you are not satisfied with the strength of Auto-Phone, you can reposition the magnet or add additional magnets.

Prior to placing the magnet on the telephone or cell phone, use a recommended cleaning agent to clean the telephone.

How to use Auto-Phone

1. Lift the telephone to your ear.
2. When you hear a short melody, the phone program is active.

NOTE:

- You may need to move the telephone receiver slightly to find the best position for a reliable Auto-Phone activation and a good hearing experience on the telephone.
- If your hearing aids have enabled the Asymmetric Phone Handling functionality, the hearing aid on the non-phone ear automatically attenuates.
- Do not cover the phone loudspeaker opening with the magnet.
- If the function does not work to your satisfaction, moving the magnet to another position may improve ease of use and comfort.
- If your hearing aids do not switch to the telephone program consistently, try repositioning the magnet or adding additional magnets.
- Use a recommended cleaning agent.

Advanced options

Using your hearing aids with iPhone, iPad, and iPod touch (optional)

The advanced models of our hearing aids are Made for iPhone, iPad and iPod touch, which allow for direct audio streaming and control from these devices.

Beltone Remote Care and Beltone Remote Care Live (optional)

Beltone Remote Care

If you have signed up to use Beltone Remote Care available with your hearing aids, you can allow your hearing aids to be adjusted remotely without having to visit your hearing care professional.

All you need is a compatible mobile device with internet enabled. This allows you to:

- Request assistance remotely to adjust your hearing aids to be a better fit for you.
- Keep your hearing aids up to date with the latest software to ensure the best performance possible.

This service only works if your mobile device is connected to the internet. Your hearing care professional will provide information regarding this option, and how it works with the Beltone HearMax™ app.

For optimum performance, make sure the hearing aids are connected to the Beltone HearMax™ app and placed close to the iPhone, iPad, iPod touch, or the Android™ smartphone before applying the changes.

NOTE: Your hearing aids shut down during the installation and update process.

Beltone Remote Care Live

This service also includes Beltone Remote Care Live. With this service you can get face-to-face assistance from your hearing care professional from home.

Streaming from an Android™ smartphone

Some Android smartphones can stream audio directly to the advanced models of our hearing aids. Your device must be running Android 10 or newer and it must have the Android Streaming for Hearing Aids feature as well.

NOTE: For assistance with pairing and using these products with your hearing aids, contact your hearing care professional.

Using your hearing aids with smartphone apps (optional)

The smartphone apps send and receive signals to and from the hearing aids via smartphones.

- Do not disable app notifications.
- Install updates to keep the app working correctly.
- The app must only be used with Beltone hearing aids for which they are intended, and Beltone takes no responsibility if the app is used with other hearing aids.
- If you want a printed version of the user guide for the app, please go to our website (see the back page of this user guide) or consult customer support.

NOTE:

- For assistance with pairing and using these products with your hearing aids, contact your hearing care professional or visit our support site.
- If your Bluetooth® enabled Android smartphone does not stream directly to your hearing aids, you can use our Beltone Direct Phone Link 2 for streaming capabilities.



Flight Mode (optional)

Your hearing aids can be controlled from your mobile device or Remote Control, an option that can be added by your hearing care professional. However, in some areas you are requested to turn off wireless communication.



CAUTION: When boarding a flight or entering an area where RF transmitters are prohibited, wireless functionality must be deactivated.

NOTE: Both hearing aids must be set in Flight mode - even with synchronisation enabled.

Turning off wireless communication (enter Flight Mode)

1. Open and close the battery door on each hearing aid three times within 10 seconds.
2. A 10-second double tone (♪♪) means the hearing aid is now in Flight Mode.

NOTE: Both hearing aids must be set in Flight mode - even with synchronisation enabled.

Activating wireless communication (exit Flight Mode)

1. Open and close the battery door on each hearing aid once.
2. Wireless communication will be activated after 10 seconds.

NOTE: It is important to wait an additional 15 seconds after wireless function resumes before opening and closing the battery compartment again for any reason. Flight mode will resume if you open and close the battery compartment during this 15 second window.

Wireless accessories

Our wireless eco-system features a comprehensive range of seamlessly integrated wireless accessories. These let you stream high quality stereo sound and speech directly to your hearing aids.

- **A TV streamer** streams audio from a TV and most other audio sources to your hearing aids at a volume that suits you.
- **A basic remote control** adjusts volume, mutes your hearing aids and changes programs.
- **An advanced remote control** adjusts volume, mutes your hearing aids, changes programs and displays your settings.
- **A phone clip** streams phone conversations and stereo sound to both hearing aids and doubles as a remote control.
- **A body-worn microphone** is a microphone that can be worn by others. It improves speech comprehension in noisy situations.
- **A wireless microphone.** It works like the **body-worn microphone**, but doubles as a table microphone. Furthermore, it has a built-in telecoil that allows it to connect with a teleloop system, a connector for a FM receiver, and a mini-jack input for wired streaming of audio from a computer or music player.

Accessing wireless accessories

To access a wireless accessory that has already been paired with your hearing aids, press the push button for 2 seconds. The hearing aid will emit a sound to confirm the connection.

For information about how to pair your hearing aids with a wireless accessory, see the user guide for the relevant wireless accessory.

NOTE:

- Ask your hearing care professional for more information on the range of Beltone wireless accessories.
- You should only use Beltone wireless accessories with your wireless hearing aids.

How to clean and maintain your hearing aids

Daily maintenance

NOTE: It is important to keep your hearing aid clean and dry. On a daily basis, clean the hearing aids using a soft cloth or tissue. In order to avoid damage due to humidity or excessive perspiration, the use of a drying kit is recommended.

NOTE: Do not use alcohol or any other solvents to clean your hearing aids - the protective coating will be damaged.

Care and maintenance

Please follow the advice below to have the best user experience and to prolong the life of your hearing aids.

- Keep your hearing aids dry and clean.
- Open the battery door to dry your hearing aids when you are not wearing them.
- Wipe the hearing aids with a soft cloth after use to remove grease or moisture.
- Do not wear your hearing aids when putting on cosmetics, perfume, aftershave, hairspray, suntan lotion, etc. These might discolour the hearing aid or get into the hearing aid causing damage.
- Do not immerse your hearing aid in any liquid.
- Keep your hearing aids away from excessive heat and intense direct sunlight. The heat may deform the shell, damage the electronics and deteriorate the surface.
- Do not swim, shower or steam bathe while wearing your hearing aids.

How to replace wax filters

Custom hearing aids may have wax filters that protect against wax and moisture. It is recommended that these are changed as needed.

Change HF3 wax filters

1. Brush the sound outlet area with the sound outlet pointed down.
2. Insert the threaded end of the wax filter tool into the used wax filter, and gently rotate clockwise.
3. Gently pull until the used filter is removed.
4. Discard the used filter in the slot located in the wax filter kit by pressing it into the centre, sliding it to one end of the slot, and pull until the filter is discarded.
5. Flip the wax filter tool around, locate a new filter in the dial, and press the tip of the tool into the centre of the dial.
6. Gently pull the new filter out of the dial.
7. Align the new filter to the sound outlet.
8. Press the new filter into the opening, and simultaneously pull and rock back and forth until the new wax filter is in place.

Change wax filters

1. To remove the old wax guard, insert the removal end of the wax guard tool into the used wax guard so that the shaft of the tool is touching the rim of the wax guard. Slowly pull the wax guard straight out.
2. To insert the new wax guard, gently press the replacement end of the wax guard tool straight into the hole of the sound outlet until the outer ring lies flush with the outside of the receiver. Pull the tool straight out - the new wax guard should remain in place.

NOTE: Pressing on the new filter with the flat side of the wax filter tool can ensure that the filter is correctly in place.

If a different type of wax filter is used for your hearing aids, or if your hearing aids do not use wax filters, consult your hearing care professional for guidance.



CAUTION: Use only original consumables from the manufacturer e.g. wax filters.

General warnings and precautions



General warnings

1. Consult a hearing care professional if you think there may be a foreign object in your ear canal, if you experience skin irritation, or if excessive earwax accumulates with the use of the hearing aid.
2. Different types of radiation, from e.g. NMR, MRI, or CT scanners, may damage hearing aids. It is recommended not to wear hearing aids during these or other similar procedures. Other types of radiation, such as burglar alarms, room surveillance systems, radio equipment, mobile telephones contain less energy and will not damage hearing aids. However, they have the potential to momentarily affect the sound quality or temporarily create undesired sounds from the hearing aids.
3. Do not wear hearing aids in mines, oil fields, or other explosive areas unless those areas are certified for hearing aid use.
4. Do not allow others to use your hearing aids.
5. Hearing aid usage by children or mentally disabled persons should be supervised at all times to ensure their safety. The hearing aid contains small parts that could be swallowed by children. Please be mindful not to leave children unsupervised with this hearing aid.
6. Hearing aids should be used only as prescribed by your hearing care professional. Incorrect use may result in sudden and permanent hearing loss.
7. Warning to hearing care professionals: Special care should be exercised in selecting and fitting hearing aids with maximum sound pressure level that exceeds 132 dB SPL with an IEC 60711:1981 occluded ear simulator. There may be a risk of impairment of the remaining hearing.

8. Turn off your wireless functionality by using the flight mode in areas where radio frequency emission is prohibited.
9. If a hearing aid is broken, do not use it.
10. A power hearing aid can produce very loud sound to compensate for a severe or profound hearing loss. There is therefore a risk of further impairing the remaining hearing.
11. External devices connected to the electrical input must be safe according to the requirements of IEC 60601-1, IEC 60065, EN/IEC 62368-1, or IEC 60950-1, as appropriate (wired connection, for example HI-PRO, SpeedLink).

NOTE:

For use of wireless functionality, only use supported wireless accessories. For further guidance regarding pairing etc., please refer to the user guide of the relevant wireless accessory.



General precautions

- When wireless function is activated, the device uses low-powered digitally coded transmissions in order to communicate with other wireless devices. Although unlikely, nearby electronic devices may be affected. In that case, move the hearing aid away from the affected electronic device
- Use only original parts from the manufacturer, e.g. wax guards.
- Only connect your hearing aids to accessories intended and qualified to be used with your hearing aids.

Hearing aid expectations

- A hearing aid will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions.
- Consistent use of the hearing aid is recommended. In most cases, infrequent use does not permit you to get full benefit from it.
- The use of a hearing aid is only part of hearing rehabilitation and may need to be supplemented by auditory training and instructions in lip-reading.

Troubleshooting

Issue	Potential cause	Potential solution
Feedback, "whistling"	Is the hearing aid inserted correctly?	Re-insert it.
	Is the volume very loud?	Reduce it.
	Are you holding an object (e.g. a hat, a telephone receiver) close to the hearing aid?	Move your hand away to create more space between the hearing aid and the object.
	Is your ear full of wax?	Visit your physician.
No sound	Is the hearing aid turned on?	Turn it on.
	Is there a battery in the hearing aid?	Insert a new battery.
	Is the battery still good?	Replace with a new one.
	Is your ear full of wax?	Visit your physician.
Sound is distorted or weak?	Did your hearing aid get moist?	Use a desiccant (drying kit).
	Is the battery dirty?	Clean it or replace it with a new one.
	The battery is dead	Replace it with a new one.
Battery drains very quickly	Did you leave your hearing aid turned on for long periods of time?	Always switch off your hearing aids when you are not using them, e.g. during the night.
	Is the battery old?	Check the date on the battery pack.

Tinnitus Management

Tinnitus Breaker Pro module

Your Beltone hearing aids include the Tinnitus Breaker Pro (TBP) module. The Tinnitus Breaker Pro (TBP) Module is a software tool that generates sounds to be used in tinnitus management programs to relieve suffering from tinnitus. The TBP can generate sounds adjusted to your personal preference and your specific therapeutic needs as determined by your doctor, audiologist or hearing care professional. Depending on the selected hearing aid program and the environment you are in, you may hear the therapeutic sound of a continuous or fluctuating noise.

User instructions for the TBP module

Description of the device

The Tinnitus Breaker Pro (TBP) Module is a software tool that generates sounds to be used in tinnitus management programs to temporarily relieve suffering from tinnitus.

Explanation of how the device works

The TBP module is a frequency and amplitude shaped white-noise generator. The noise signal level and frequency characteristics can be adjusted to the specific therapeutic needs as determined by your doctor, audiologist or hearing care professional.

Your doctor, audiologist or hearing care professional can modulate the generated noise with the purpose of making it more pleasant. The noise can then resemble, for example, breaking waves on a shore.

Modulation level and speed can also be configured to your likes and needs. An additional feature can be enabled by your hearing care professional that allows you to select predefined sounds that simulate sounds from nature, such as breaking waves or running water.

If you have two wireless hearing aids that support ear-to-ear synchronisation, this functionality can be enabled by your hearing care professional. This will cause the Tinnitus Breaker Pro to synchronise the sound in both hearing aids.

If your tinnitus troubles you only in quiet environments, your doctor, audiologist or hearing care professional can set the TBP Module so that it becomes audible exclusively in such surroundings. The overall sound level can be adjusted via a volume control. Your doctor, audiologist or hearing care professional will review with you the need for having such a control.

For hearing aids where ear-to-ear synchronisation is enabled, your hearing care professional can also enable environmental monitoring synchronisation so that the TBP noise level is automatically adjusted simultaneously in both hearing aids dependent on the background sound level. Additionally, since the hearing aid has a volume control, the background noise level monitored by the hearing aid and the volume control can be used simultaneously to adjust the generated noise level in both hearing aids.

The scientific concepts that form the basis for the device

The TBP module provides sound enrichment with the aim of surrounding the tinnitus sound with a neutral sound which is easily ignored. Sound enrichment is an important component of most approaches to tinnitus management, such as Tinnitus Retraining Therapy (TRT).

To assist habituation to tinnitus, this needs to be audible. The ideal level of the TBP module, therefore, should be set so that it starts to blend with the tinnitus, and so that you can hear both your tinnitus as well as the sound used.

In a majority of instances, the TBP module can also be set to mask the tinnitus sound, so to provide temporary relief by introducing a more pleasant and controllable sound source.

TBP volume control

The sound generator is set to a specific loudness level by the hearing care professional. When switching the sound generator on, the volume will have this optimal setting. Therefore, it might not be necessary to control the volume (loudness) manually. However, the volume control provides the ability to adjust the volume, or amount of stimulus, to the liking of the user. The tinnitus sound generator volume can only be adjusted within the range set by the hearing care professional.

The volume control is an optional feature in the TBP module used for adjusting the sound generator output level.

Using TBP with smartphone apps

The tinnitus sound generator control via hearing aid push buttons can be enhanced with wireless control from a TBP control app on a smartphone or mobile device. This functionality is available in supported hearing aids when a hearing care professional has enabled the TBP functionality during fitting of the hearing aid.

NOTE: To use smartphone apps, the hearing aid must be connected with the smartphone or mobile device.

TBP - Technical specifications

Audio signal technology

Digital.

Available sounds

White noise signal which can be shaped with the following configurations:

High-pass filter	Low-pass filter
500 Hz	2000 Hz
750 Hz	3000 Hz
1000 Hz	4000 Hz
1500 Hz	5000 Hz
2000 Hz	6000 Hz
-	8000 Hz

The white noise signal can be modulated in amplitude with an attenuation depth of up to 14 dB.



Prescription use of a Tinnitus Breaker Pro hearing aid

The TBP should be used as prescribed by your doctor, audiologist or hearing healthcare professional. In order to avoid permanent hearing damage, the maximum daily usage depends on the level of the generated sound.

To adjust TBP, please consult your hearing care professional.

Should you develop any side effects from using the sound generator, such as dizziness, nausea, headaches, perceived decrease in auditory function or increase in tinnitus perception, you should discontinue use of the sound generator and seek medical evaluation.

Children and physically or mentally challenged users will require training by a doctor, audiologist, hearing care professional or the guardian for the insertion and removal of the hearing aid containing the TBP module.

Important notice for prospective sound generator users

A tinnitus masker is an electronic device intended to generate noise of sufficient intensity and bandwidth to mask internal noises. It is also used as an aid in hearing external noises and speech.

Good health practice requires that a person with a tinnitus condition have a medical evaluation by a licensed physician (preferably a physician who specializes in diseases of the ear) before using a sound generator. Licensed physicians who specialize in diseases of the ear are often referred to as otolaryngologists, otologists or otorhinolaryngologists.

The purpose of medical evaluation is to assure that all medically treatable conditions that may affect tinnitus are identified and treated before the sound generator instrument is used.

The sound generator instrument is a tool to generate sounds to be used with appropriate counselling and/or in a tinnitus management program to relieve patients suffering from tinnitus.

Warning information



WARNING:

- Sound generators can be dangerous if improperly used.
- Sound generators should be used only as advised by your doctor, audiologist, or hearing care professional.

- Sound generators are not toys and should be kept out of reach of anyone who might cause themselves injury (especially children and pets).



CAUTION:

- Should the user develop any side effects from using the sound generator, such as dizziness, nausea, headaches, perceived decrease in auditory function or increase in tinnitus perception, the user should discontinue use of the sound generator and seek medical evaluation.
- To prevent unintended usage by paediatric or physically or mentally disabled users, the volume control must, if enabled, be configured to only provide a decrease of the sound generator output level.
- Children and physically or mentally disabled users will require guardian supervision while wearing the TBP hearing aid.



Tinnitus Breaker Pro warning to hearing care professionals

A hearing care professional should advise a prospective sound generator user to consult promptly with a licensed physician (preferably an ear specialist) before getting a sound generator, if the hearing care professional determines through inquiry, actual observation, or review of any other available information concerning the prospective user that the prospective user has any of the following conditions:

1. Visible congenital or traumatic deformity of the ear.
2. History of active drainage from the ear within the previous 90 days.
3. History of sudden or rapidly progressive hearing loss within the previous 90 days.
4. Acute or chronic dizziness.

5. Unilateral hearing loss of sudden or recent onset within the previous 90 days.
6. Audiometric air-bone gap equal to or greater than 15 dB at 500 Hertz (Hz), 1000 Hz, and 2000 Hz.
7. Visible evidence of significant cerumen accumulation or a foreign body in the ear canal.
8. Pain or discomfort in the ear.



CAUTION: The maximum output of the sound generator falls into the range that can cause hearing loss according to OSHA regulations. In accordance with NIOSH recommendations, the user should not use the sound generator for more than eight (8) hours a day when this is set to a level of 85 dB SPL or above. When the sound generator is set to levels of 90 dB SPL or above the user should not use the sound generator for more than two (2) hours per day. In no case should the sound generator be worn at uncomfortable levels.



Tinnitus Breaker Pro precautions

1. Should the user develop any side effects from using the sound generator, such as dizziness, nausea, headaches, perceived decrease in auditory function or increase in tinnitus perception, the user should discontinue use of the sound generator and seek medical evaluation.
2. Discontinue use of the sound generator and consult promptly with a licensed physician if you experience one of the following conditions:
 - a. Visible congenital or traumatic deformity of the ear.
 - b. History of active drainage from the ear within the previous 90 days.
 - c. History of sudden or rapidly progressive hearing loss within the previous 90 days.
 - d. Acute or chronic dizziness.
 - e. Unilateral hearing loss of sudden or recent onset within the previous 90 days.

- f. Visible evidence of significant cerumen accumulation or a foreign body in the ear canal.
 - g. Pain or discomfort in the ear.
3. Discontinue use of the sound generator and consult promptly with your hearing care professional, if you experience changes in the tinnitus perception, discomfort or interrupted speech perception, while using the Tinnitus Breaker Pro.
 4. The volume control is a feature in the TBP module used for adjusting the sound generator output level. To prevent unintended usage by paediatric or physically or mentally disabled users, the volume control must be configured to only provide a decrease of the sound generator output level.
 5. Children and physically or mentally disabled users will require guardian supervision while wearing the TBP hearing aid.
 6. Adjustment of the Tinnitus Breaker Pro settings, using a smartphone app, should only be performed by the parent or legal guardian in cases where the user is minor. Use of the [[[Undefined variable CrossBrand_FeatureNames.epre]]] for remote settings of the tinnitus sound generator, should only be performed by the parent or legal guardian in cases where the user is minor.

Regulatory information

Warranties and repairs

The manufacturer provides a warranty on hearing aids in the event of defects in workmanship or materials, as described in applicable warranty documentation. In its service policy, the manufacturer pledges to secure functionality at least equivalent to the original hearing aid. As a signatory to the United Nations Global Compact initiative, the manufacturer is committed to doing this in line with environment-friendly best practices. Hearing aids therefore, at the manufacturer's discretion, may be replaced by new products or products manufactured from new or serviceable used parts, or repaired using new or refurbished replacement parts. The warranty period of hearing aids is designated on your warranty card, which is provided by your hearing care professional.

For hearing aids that require service, please contact your hearing care professional for assistance.

Hearing aids that malfunction must be repaired by a qualified technician. Do not attempt to open the case of hearing aids, as this will invalidate the warranty.

Temperature test, transport and storage information

Our hearing aids are subjected to various tests in temperature and damp heating cycling between -25 °C (-13 °F) and +70 °C (+158 °F) according to internal and industry standards.

During normal operation the temperature should not exceed the limit values of 0 °C (+32 °F) to +45 °C (+113 °F) at a relative humidity of 90%, non-condensing. An atmospheric pressure between 500 hPa and 1100 hPa is appropriate.

NOTE: During use, your hearing aids may reach temperatures up to +43 °C (+109 °F).

During transport or storage, the temperature should not exceed the limit values of -20 °C (-4 °F) to +60 °C (+140 °F) at a relative humidity of 90% RH, non-condensing (for a limited time).

Statement

This device complies with part 15 of the FCC rules and ISSED rules. Its operation is subject to the following two conditions:

1. This device may not cause harmful interference.
2. This device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC rules and ISSED rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used according to the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment and the receiver

- Connect the equipment to an outlet or a circuit that is different from the one to which the receiver is connected
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications can void the user's authority to operate the equipment.

The products are in compliance with the following regulatory requirements:

- In the EU: The device conforms to the Essential Requirements according to Annex I of Council Directive 93/42/EEC for medical devices (MDD).
- Hereby, Beltone A/S declares that the radio equipment types CSX10, CSI12 and CSI13 are in compliance with Directive 2014/53/EU.
- The full text of the EU declaration of conformity is available at the following internet address: www.beltone.com/en/declarations.
- Other identified applicable international regulatory requirements in countries outside the EU and US. Please refer to local country requirements for these areas.
- In Canada: these hearing aids are certified under the rules of ISED.
- Japanese Radio Law and Japanese Telecommunications Business Law Compliance: This device has been certified pursuant to the Japanese Radio Law (電波法) and the Japanese Telecommunications Business Law (電気通信事業法). This device should not be modified (otherwise the granted designation number will become invalid).

Type designations

Hearing aid type designations for models included in this user guide are:

CSI12, FCC ID: X26CSI12, IC: 6941C-CSI12 **CSI13**, FCC ID: X26CSI13, IC: 6941C-CSI13, and **CSX10**, FCC ID: X26CSX10, IC: 6941C-CSX10.

This device includes an RF transmitter which operates in the frequency band of 2.4 GHz – 2.48 GHz.

Hearing aid variantsSome models may not be available in your region.

Completely-in-the-canal (CIC) hearing aids (including type CSX10 with FCC ID X26CSX10, IC number 6941C-CSX10 models) with size 10A battery, are available in the following variants:

RLY4CIC-W-LP, RLY3CIC-W-LP,
RLY4CIC-W-MP, RLY3CIC-W-MP,
RLY4CIC-W-HP, RLY3CIC-W-HP

In-the-canal (ITC) hearing aids (including type CSI12 with FCC ID X26CSI12, IC number 6941C-CSI12 models) with size 312 battery, are available in the following variants:

RLY4ITC-DW-MP, RLY3ITC-DW-MP, RLY2ITC-DW-MP,
RLY4ITC-DW-HP, RLY3ITC-DW-HP, RLY2ITC-DW-HP,
RLY4ITC-DW-UP, RLY3ITC-DW-UP, RLY2ITC-DW-UP

In-the-ear (ITE) hearing aids (including type CSI13 with FCC ID X26CSI13, IC number 6941C-CSI13 models) with size 312 or 13 battery, are available in the following variants:

RLY4ITE-DW-MP, RLY3ITE-DW-MP, RLY2ITE-DW-MP,
RLY4ITE-DW-HP, RLY3ITE-DW-HP, RLY2ITE-DW-HP,
RLY4ITE-DW-UP, RLY3ITE-DW-UP, RLY2ITE-DW-UP

Nominal RF output power transmitted is :

CSI12: + 1.36 dBm - **CSI13:** + 1.36 dBm - **CSX10:** + 2.5 dBm.

Symbols



WARNING: Points out a situation that could lead to serious injuries.



Legal manufacturer.



CAUTION: Indicates a situation that could lead to minor and moderate injuries.



Date of manufacture.



Follow instructions for use. (Logo in blue)



Product is a Type B applied part.



Unique Device Identification.



Equipment includes an RF transmitter.



IMDA label for Singapore.



Regulatory compliance mark for Australia and New Zealand.



Do not dispose of electronic devices or batteries with ordinary household waste. They should be disposed of at sites intended for electronic waste or returned to your hearing care professional for safe disposal.

Please note: There may be specific regulations in your country.

Additional information

Technical specifications

Hearing aid model	Maximum output (Data in accordance with IEC60118-0 Edition 3.0 2015-06, IEC60118-7 and ANSI S3.22-2009, supply Voltage 1.3V)
All Low Power (LP) models (CIC models)	114 dB SPL (typical)
All Medium Power (MP) models	118 dB SPL (typical)
All High Power (HP) models	120 dB SPL (typical)
All Ultra Power (UP) models (Not CIC models)	129 dB SPL (typical)

Acknowledgements

Portions of this software are written by Kenneth MacKay (micro-ecc) and licensed under the following terms and conditions:

Copyright © 2014, Kenneth MacKay. All rights reserved.

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

- Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
- Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.

THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT HOLDER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

NOTE: Use of the Made for Apple badge means that an accessory has been designed to connect specifically to iPhone, iPad, and iPod touch models, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards.

© 2024 GN Hearing A/S. All rights reserved. Beltone is a trademark of GN Hearing A/S. Apple, the Apple logo, iPhone, iPad, iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the US and other countries. Android, Google Play and the Google Play logo are trademarks of Google LLC. The Bluetooth word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

Notes

Notes

Notes

Manufacturer according to EU Medical
Device Directive 93/42/EEC:



Beltoné A/S
Lautrupbjerg 7
DK-2750 Ballerup
Denmark
beltone.com

United Kingdom

GN Hearing UK Ltd
Unit 13 Talisman Business Centre
Bicester OX26 6HR
United Kingdom
Tel.: +44 1869 352 800
beltone.com

Singapore

GN Hearing Pte. Ltd
456 Alexandra Road
Fragrance Empire Building #22-01
Singapore 119962
Tel: +65 6320 9388
beltone.com

Australia

GN Hearing Australia Pty Ltd
Gate C, 19-25 Khartoum Rd
Macquarie Technology Park
Macquarie Park NSW 2113
Australia
Tel.: (free) 1800 658 955
beltone.com

New Zealand

GN Hearing New Zealand Limited
Ground Floor, North Entrance
4 Fred Thomas Drive
Takapuna, Auckland, 0622
Tel.: (free) 0800 900 126
beltone.com



Any serious incident that has occurred in relation to the device should be reported to the Legal manufacturer Beltoné A/S and the competent authority of the EU Member State in which the user and/or patient is established.